

HIGH PERFORMANCE EMCCD & CCD CAMERAS FOR LIFE SCIENCES



Primary applications
Live-cell imaging
High-speed emission ratio imaging
Low-copy gene analysis and gene expression profiling
Quantitative FRET, FRAP, FISH
Luminescence






1392 x 1040 imaging array 6.45 x 6.45-µm pixels

The CoolSNAP HQ2 Monochrome camera from Photometrics® delivers fast, high-resolution imaging for quantitative fluorescence microscopy applications. This cooled CCD camera provides a large dynamic range with very low noise at both 10 MHz and 20 MHz. The fine pitch of the pixels is ideally matched to the resolution of optical microscopes. Megapixel resolution and small pixels allow imaging of very fine detail, yet the pixels can be easily binned to improve sensitivity. Advanced interline-transfer CCD technology provides high quantum efficiency, most notably in the near-infrared (NIR) portion of the spectrum.

Features	Benefits
10-MHz and 20-MHz readout	Dual-mode readout for variable-speed image capture
1392 x 1040 imaging array 6.45 x 6.45-µm pixels	Resolves fine detail Ideally matched to optical microscope
Interline-transfer, progressive-scan CCD	Electronic shuttering eliminates camera vibration and facilitates fast triggering
Flexible binning and readout	Increases signal-to-noise performance while increasing the frame rate
IEEE-1394a	High-bandwidth, uninterrupted data transfer with no dropped frames
Digitization	Quantifies bright and dim signals in the same image 14-bit digitization
Thermoelectric cooling	Special cooling package virtually eliminates dark current
Enhanced quantum efficiency	Provides higher sensitivity than typical interline cameras (especially in the NIR)
C-mount	Easily attaches to microscopes, standard lenses, or optical equipment
Acquisition software	Captures, analyzes, and saves high-resolution images
PVCam® Circular buffers Device sequencing	Supported by numerous third-party software packages Real-time focus Precise integration with shutters, filter wheels, etc.
IEEE-1394a compatibility	Windows® XP/Vista 32

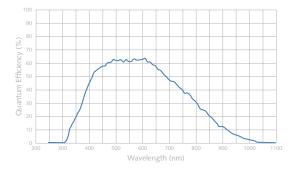

		Region			
		1392 x 1040	512 x 512	256 x 256	
Binning	1 x 1	11	21	36	
	2 x 2	20	36	58	
	3 x 3	28	48	71	
	4 x 4	35	57	81	
	8 x 8	56	81	104	
	(Frames per second)				

Note: Frame rates are measured at 20 MHz with up to 90-millisecond exposure times.

	Specifications
CCD image sensor	Sony® ICX285; interline-transfer, progressive-scan device with microlenses
CCD format	1392 x 1040 imaging array 6.45 x 6.45-µm pixels 8.98 x 6.71-mm imaging area (optically centered)
Grade	Sony Grade 0
System gain	1 e-/ADU
Linear full well	16,000 e- (single pixel) 30,000 e- (2 x 2 binned pixel)
Read noise	4.5 e- rms @ 10 MHz 5.5 e- rms @ 20 MHz
Nonlinearity	<1%
Digitization	14 bits @ 20 MHz or 10 MHz (software selectable)
Interface	IEEE 1394a FireWire
Frame readout	90 ms/frame
CCD temperature	-30°C (regulated)
Dark current	0.001 e-/p/s @ -30°C
Operating environment	0 to 30°C ambient, 0 to 80% relative humidity noncondensing
Dimensions	4.5" x 4.0" x 7.0" (6.5 lbs)
1/0	TTL (trigger/status): trigger, invert, inhibit, exposing, interline shift, frame readout 8-bit TTL (programmable) 8-bit DACs (two)

Note: Specifications are typical and subject to change.

CoolSNAP is a trademark of Photometrics. Photometrics and PVCam are registered trademarks of Photometrics. Linux is a registered trademark of Linus Torvalds. Mac OS is a trademark of Apple Computer, Inc., registered in the U.S. and other countries. Sony is a registered trademark of Sony Corporation. Windows is a registered trademark of Microsoft Corporation in the United States and other countries. Other brand and product names are the trademarks or registered trademarks of their respective owners and manufacturers.





info@photometrics.com tel: +1 520.889.9933